



KEY FEATURES

ACKERMANN XPRESS LOCK® SURGERY'S FUTURE



ACKERMANN XPRESS LOCK® DEMOUNTABLE INSTRUMENT LINE

- ✓ Latest instrument system with **quick-locking mechanism**
- ✓ Innovative ratchet design for on-off-action
- ✓ Autoclavable
- ✓ Lengths of 260, 330 and 450 mm
- ✓ 360° rotatable with audible and tactile feedback
- ✓ Demountable into handle, sheath and insert
- ✓ Most ergonomic handles in different designs for every surgeon's physique
- ✓ No hidden spaces to accommodate blood and debris
- ✓ Time saving due to easy cleaning and assembly
- ✓ Cost effective - design reduces cost per case and inventory management
- ✓ Virtually no stock if compared to disposable instruments

1

CONCEPT

The modern, detachable instrument line from Ackermann is equipped with an unique quick locking mechanism, which makes the instruments' assembly and disassembly even more efficient and faster.

2

HOW IT WORKS



ASSEMBLY AND DISASSEMBLY WITHIN 5 SECONDS !

- 1** Open the handle completely.
- 2** Loosen the sheath by pressing the quick locking knob on the black screw connector.
- 3** While holding the sheath, rotate the insert clockwise 45° until loosened.
Slide the insert from sheath.

Ackermann®

Rev. 2017-10



Ackermann Instrumente GmbH
Eisenbahnstrasse 65-67
78604 Rietheim-Weilheim
Germany



info@ackermanninstrumente.de



+49 (0)7461 966 17 - 0

www.ackermanninstrumente.de



TOPSELLERS

ACKERMANN XPRESS LOCK[®] SURGERY'S FUTURE

1

ORDERING INFORMATION



13-1461XPI
Metzenbaum,
curved,
Ø 5 mm, length 330mm

+ 13-1356XP



+ 13-1367XP



13-1431XPI
Hook scissors,
Ø 5 mm, length 330mm

+ 13-1356XP



+ 13-1367XP



13-1318XPI
Maryland,
Ø 5 mm, length 330mm

+ 13-1356XP



+ 13-1367XP



13-1323XPI
Endo Clinch,
extremely atraumatic serrated,
Ø 5 mm, length 330mm

+ 13-1356XP



+ 13-1368XP



13-1354XPI
Standard grasper,
cross serrated,
Ø 5 mm, length 330mm

+ 13-1356XP




+ 13-1368XP



13-1570XPI
De Bakey,
Ø 5 mm, length 330mm

+ 13-1356XP

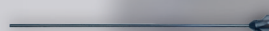

+ 13-1368XP



13-1312XPI
Johan grasping forceps,
slightly flat jaws,
Ø 5 mm, length 330mm

+ 13-1356XP




+ 13-1368XP



13-1541XPI
Babcock,
Ø 5 mm, length 330mm

+ 13-1356XP

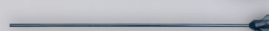

+ 13-1368XP



13-1605XPI
Claw forceps,
traumatic jaws,
Ø 10 mm, length 330mm

+ 13-1357XP



+ 13-1368XP



13-1660XPI
Big spoon forceps,
Ø 10 mm, length 330mm

+ 13-1357XP



+ 13-1367XP



Ackermann[®]

Rev. 2017-10



Ackermann Instrumente GmbH
Eisenbahnstrasse 65-67
78604 Rietheim-Weilheim
Germany



info@ackermanninstrumente.de



+49 (0)7461 966 17 - 0

www.ackermanninstrumente.de